

Optional Accessory

Infusion Bag Hook KC-0584

Accommodates large bags up to 2,000 cc in volume.

Attaching large-volume infusion bags directly to the existing hooks on the IV stand could affect the stability of the IV stand. Use of the infusion bag hook will allow you to attach infusion bags with volumes of up to 2,000 cc to the IV stand at a more stable position close to the support.

\*Dedicated option for the Paramount Bed IV Stand KC-508/KC-508A.



Urine Collection Bag Holder KC-058

Enables quick and easy attachment of urine collection bags.

Urine collection bags can be easily attached to the IV stand support using the fixing knob.

\*Dedicated option for the Paramount Bed IV Stand KC-508/KC-508A.



List of Products

Product name	Model No.	Dimensions (cm)	Product weight
IV stand	K C - 5 0 8 A	Legs: 45.3, max. diameter × 128.7 to 196.2, total height	Approx. 4.3 kg
Oxygen cylinder holder	K C - 0 5 8 1	12.7 (W) × 23.7 (D) × 41.2 (H)	Approx. 1.5 kg
Linking arm (for a bed/stretcher)	K C - 0 5 8 2	4.6 (W) × 45.8 (D) × 10.3 (H)	Approx. 500 g
Linking arm (for a wheel chair)	K C - 0 5 8 3	4.6 (W) × 45.8 (D) × 10.3 (H) (*)	Approx. 665 g
Infusion bag hook	K C - 0 5 8 4	6.3 (W) × 6.5 (D) × 11.4 (H)	Approx. 60 g
Urine collection bag holder	K C - 0 5 8	13 (W) × 13 (D) × 5.3 (H)	Approx. 250 g

\*Not including the width, depth, and height of the holder included with the stand

IV Stand KC-508A

Simple and Easy to Use



## IV Stand KC-508A

Convenient features that ensure ease of use.

IV stand embedded with a variety of extras designed to protect and help the patient.

Anti-drop mechanism to prevent sudden dropping

The support column comes equipped with an air valve designed to prevent the support column from dropping suddenly when the height is adjusted.

New

Stoppers that can be squeezed together to release and grasped together to set in place


Easy and safe height adjustment

Squeeze the upper/lower stopper grips together to adjust the height and let them go to lock each stopper in place, for easy and safe operation.




Step that can be stepped on to secure the IV stand in place

A step is provided around the legs that can easily be depressed with a foot to keep the IV stand from moving, which is useful during infusion bag replacement.

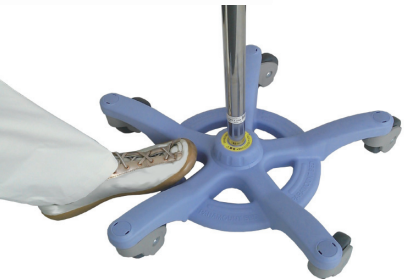

Grip that can be adjusted and grasped onto in a variety of ways


This grip, created based on universal design principles, can be adjusted as needed to accommodate the patient's condition utilizing its unique shape and angle-adjusting function.

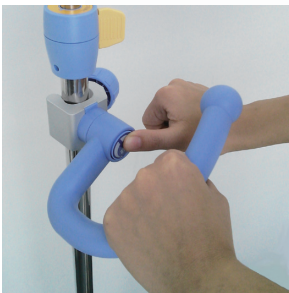

High-level traveling performance and detailed conveniences designed to benefit the patient

These castors run smoothly with little resistance even when first moved and turn on a dime. The wheel covers on each side are equipped with protrusions designed to keep the stand from tipping over should the casters slip into gaps for elevator doors, for example.

The caster wheels can be removed for cleaning to cut down on any buildup of dust, dirt, and other impurities.

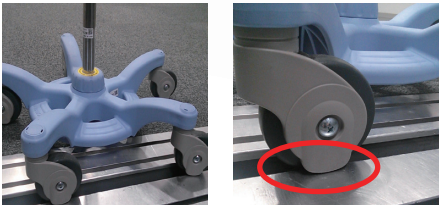

## Optional Accessory

### Oxygen Cylinder Holder KC-0581


In answer to requests from customers wanting to attach an oxygen cylinder to the IV stand

The included auxiliary stay helps to support the IV stand in order to maintain stability and keep the stand from tipping over.

It provides the same level of stability as when there is no oxygen cylinder attached to the stand, enabling you to use the IV stand safely and comfortably.

\*Dedicated option for the Paramount Bed IV Stand KC-508/KC-508A.

\*Use oxygen cylinders with a capacity of 500 L (a diameter of between 9.5 cm to 11.5 cm).

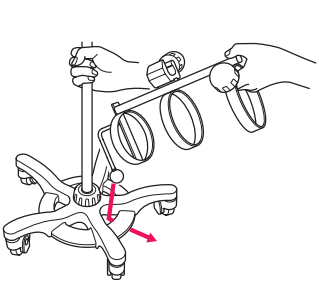

\*Insert the auxiliary stay into the hole on the leg base.




### Linking arm

Used to easily link the IV stand to other objects.


#### Bed/Stretcher Linking Arm KC-0582

The linking arm links the bed or stretcher to the IV stand at a set distance, which helps to prevent accidents caused by the IV stand moving unexpectedly and the patient's catheter getting caught on the IV stand or another object.

\*Dedicated option for Paramount Bed beds that have a hole for attaching options and for the Paramount Bed Stretcher KK-715/KK-725 Series.

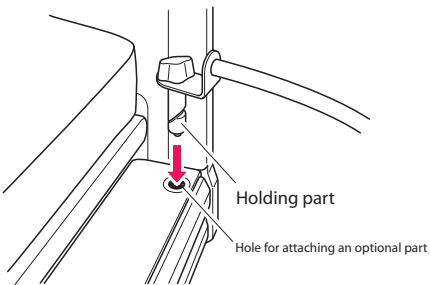

\*Use the hole for attaching options or the hole for attaching the IV hole to attach the linking arm.

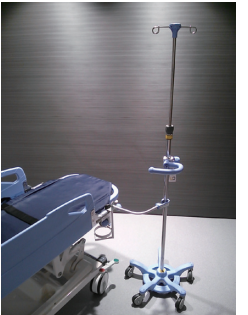

#### Wheelchair Linking Arm KC-0583

The linking arm links the wheel chair to the IV stand at a set distance allowing the IV stand to be moved safely and smoothly along side wheel chair.

As it is made of a flexible pipe material, the linking arm bends and twists to accommodate movement through narrow passageways, patient transfers from wheelchair to bed and so on, and other such situations.

\*Dedicated option for the Paramount Bed Wheelchair T50/T51/T54 Series.

\*First attach the existing holder that came with the IV stand to the wheelchair and then secure both using the linking arm.

