The materials and information provided via this brochure are distributed internationally, but not all the products or services found on this brochure are necessarily available in your country or area. Please contact local office or agents to make sure of the availability of the relevant products or services.

Mindray Medical Netherlands B.V.

Drs.W.van Royenstraat 8, 3871 AN Hoevelaken, The Netherlands Voice: (31-33) 254-4911 Fax: (31-33) 253-4280

Mindray Medical Italy S.r.l. Via Leonardo DA Vinci, 158, 20090 Trezzano Sul Naviglio (MI), Italy Voice: (39-02) 5737-401 Fax: (39-02) 5523-2018

Mindray Medical España S.L.

Avenida Manoteras 38, Edificio B, Bajo B008. 28050-Madrid Spain Voice: (34-91) 392-3754 Fax: (34-91) 392-3756

Mindray Medical Germany GmbH Zwischen den Bächen 4, 64625 Bensheim, Germany Voice: (49) 6251-17524-0 Fax: (49) 6251-17524-20

Mindray (UK) Limited

Mindray House, Kingfisher Way, Hinchingbrooke Business Park, Huntingdon, Cambs, PE29 6FN Voice: (44) 1480-416-840 Fax: (44) 1480-436-588

Mindray Medical France SARL

Europarc Creteil 1 allée des cerisiers 94000 CRETEIL, France Voice: (33-1) 4513-9150

Mindray Medical Sweden AB

Rissneleden 136, 174 57 Sundbyberg, Sweden Voice: (46) 8 555 54 121

Fax: (46) 8 555 54 101

Fax: (7-499) 553-6039

Mindray Medical Rus Co.Ltd 2 Zvenigorodskaya Street, 13, bld 41, 123022 Moscow, Russia Voice: (7-499) 553-6036

Mindray Medical Colombia S.A.S

Room No.702, Torre UNIKA, Carrera 9 No.77-67, Bogota, Colombia Voice: (57-1) 3130892 / 3210916

Mindray do Brasil - Comércio e Distribuição de Equipamentos

Médicos Ltda. Rua Tavares Bastos, 329, Perdizes, São Paulo - SP, CEP 05012-020

Voice: (55-11) 3124-8026 Fax: (55-11) 3078-8035

Mindray Medical Mexico S. de R.L. de C.V. Félix Parra # 175 Colonia. San José Insurgentes Sur Delegación Benito Juárez 03900 México, D. F. Voice: (52-55) 5661-9450 / (52-55) 5662-6620 Fax: (52-55) 5662-6597

PT.Mindray Medical Indonesia

The East 11th Floor, Unit 08 Jl. Dr. Ide Anak Agung Gde Agung Kav. E 3.2 No.1 Jakarta 12950 Voice: (6221) 2902-7280 Fax: (6221) 2902-7283

Mindray Medical Egypt Ground Floor - 11 El Imam Metwally El Sharawi Street, Sheraton Heliopolis - El Nozha, Cairo, Egypt Voice: (2) 02 22691047 Fax: (2) 02 22691046

Mindray Medical Australia

ADD: 19 Harker St. Burwood VIC 3125, Australia

Mindray Medical Nigeria

ADD: 1004 Estate, Victoria Island, Lagos, Nigeria Voice: (+234) 7030682053

Mindray Medical India Pvt. Ltd. (Delhi Office)

16th floor, Building No. 9-B, DLF Cyber City, DLF Phase- 3, Gurgaon- 122002, Haryana Voice: (91-124) 4632-488 Fax: (91-124) 4632499

Mindray Medical India Pvt. Ltd. (Mumbai office) B-404, City Point, Andheri Kurla Road, Andheri East, Mumbai, 400 059, India Voice: (91-22) 4020-0000

Mindray Medikal Teknoloji Istanbul Ltd. Sti.

KÖYALTI MEVKİİ 29 EKİN GƏN. 1 KUYUMCUKENT 2.NO.PLZ.3.KT.
D.1-9 34520 BAHÇELİEVLER/YENİBOSNA, İSTANBUL, TÜRKİYE
Voice: (90-212) 482 08 77 Fax: (90-212) 482 08 78

Mindray Medikal Teknoloji Istanbul Ltd. Sti. Branch Office Ankara TURAN GÜNEŞ BULVARI ALEKSANDER DÜBÇEK CAD. NO:18/A CANKAYA/ANKARA, TÜRKİYE Fax: (90-312) 441 96 93

Mindray Medical Vietnam Co., Ltd. Unit 2, 4th Floor, "B" Building, Broadway Office Park, 102 Nguyen Luong Bang Blvd., Tan Phu Ward, Dist. 7, HCMC, Vietnam Voice: (84) 8 5413-6970 Fax: (84) 8 5413-6971

Mindray Medical Thailand Limited 1768 Thai summit Tower Room 407, 4th floor, New petchaburee Road, Kwaeng Bangkapi, Khet Huai Khwang, Bangkok 10320, Thailand




Mindray Building, Keji 12th Road South, High-tech Industrial Park, Nanshan, Shenzhen 518057, P.R. China Tel: +86 755 8188 8998 Fax: +86 755 26582680 E-mail: intl-market@mindray.com www.mindray.com

Mindray is listed on the NYSE under the symbol "MR"

mindray | Nouthcare within reach are registered trademarks or trademarks owned by Shenzhen Mindray Bio-medical Electronics Co., LTD. © 2014 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. Specifications subject to changes without prior notice. Pls refer to user manual for manufacturer information P/N:ENG-HyBase 8300/8500-210X285X16P-20141028



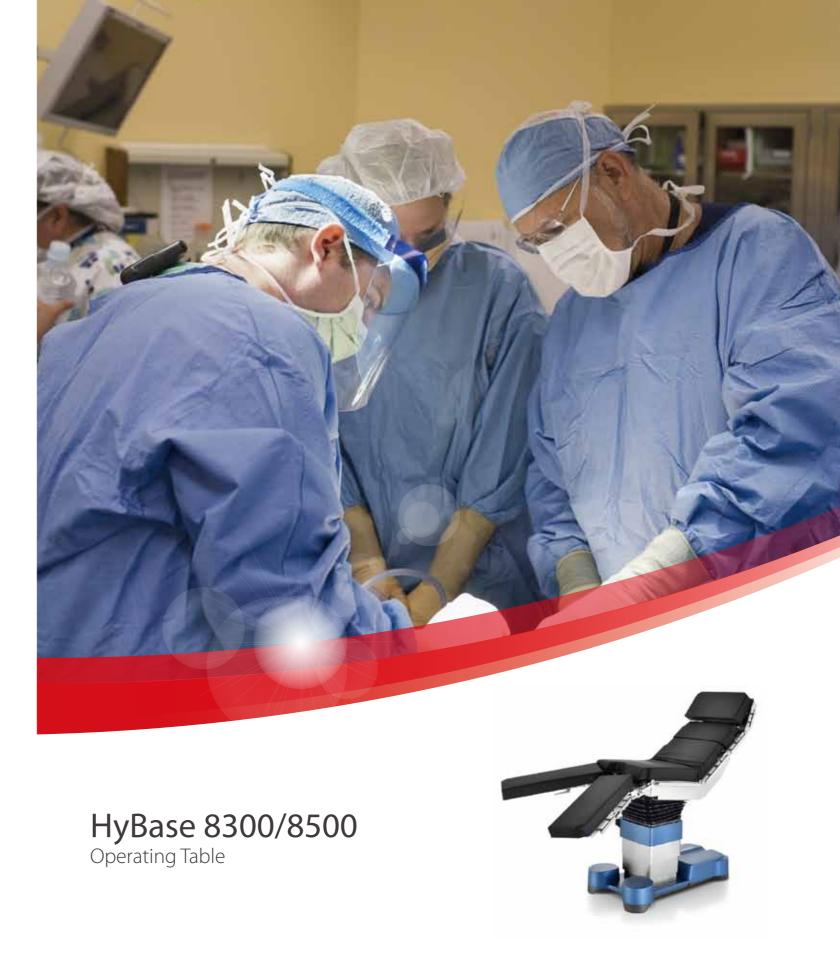







Engineered for increasingly complex surgical scenes, HyBase 8300/8500 comes with extreme weight capacity up to 460 kg guaranteeing highest safety and stability in any position. With longitudinal slide up to 480 mm, it is ideally suitable for cardiovascular surgery with C-arm. The modular design tabletop enables it to be tailored for various surgical discipline needs.

- Longitudinal slide up to 480 mm for free access to C-arm
- Extreme weight load capacity of 460 kg with highest safety and stability
- Modular design tabletop tailored for various surgical discipline needs
- Intelligent anti-collision system to avoid table components collision during movement
- Multi-layer decompression pad with antistatic, waterproof and seamless design
- 575 mm lowest table position for Neurosurgery
- IR wireless remote control system for free movements in OR

HyBase 8300/8500 Electro-Hydraulic Operating Table



Ergonomic and humanization design

- Soft TPU material coating for comfortable
- Press any key to start and auto-power off
- Backlight for operating in a dark environment
- Flex, reflex, beach chair, and "0" position with one button for each position

Better visualization - 2.4" color LCD

The useful information will be shown on the display, such as movements, battery, brake, anti-collision informa02

10 preset table positions

The table can memorize 10 positions, so the user can select the preset positions for the specific surgery freely without complex adjustment.

Intuitively UI design

The hand control button for table component movement is positioned in the corresponding patient figure chart, which enables the user to find the button fast and precisely.

Wireless charging function

Without any charging interface, the hand control is easy to clean, matching the hygiene OR standard.

05

06

IR wireless remote control system

With wireless hand control, the user can control the table in any place of the OR.

Electro-Hydraulic Operating Table



Stable brake with 4 auto-up-down castors design

Extreme Weight Capacity

HyBase 8300/8500 offers highest level of safety and stability. With a maximum weight capacity of 460 kg, the table is designed to accomodate the lightest to the heaviest surgical patient in any surgery position, ideally suitable for bariatric patients.


Intelligent Anti-collision System

In the connection joint of leg plate, there is a intelligent module that can recognize each type of the table component and match the anti-collision data automatically. So the intelligent anti-collision system can monitor the positions of the table tops during movements such as sitting plate, leg plates and so on, to avoid table components collision.



Large Longitudinal Slide

With longitudinal slide up to 480 mm, HyBase 8300/8500 enables optimum access with the C-arm without reposition the patient for imaging purpose, ideally suitable for such as cardiovascular surgery. By the way, the carbon-fibre table top (optional) provides excellent imaging quality.

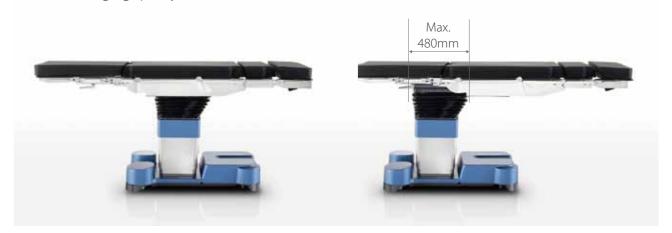

Multi-layer Decompression Pad

With the Special Foam Core pad and multi-layer design, HyBase 8300/8500 provides excellent pressure distribution, significantly avoids the main cause of decubitus, reduces the rate of postoperative health risk, and speeds up the patient recovery. Meanwhile, the pad is waterproof and easy to clean as ultrasonically-welded technology and GORE-TEX material.







Lowest Table Height

With lowest height up to 575mm for neurosurgery, the doctor can sit and do the surgery with comfort.





Electro-Hydraulic Operating Table

Modular design


Modular design table top is tailored for various surgical discipline needs, which also can low down the cost of the hospital as they only need to buy the specific module for more surgeries not one more table. With Module Recognition System, the table can recognize the module type and match movement data automatically.









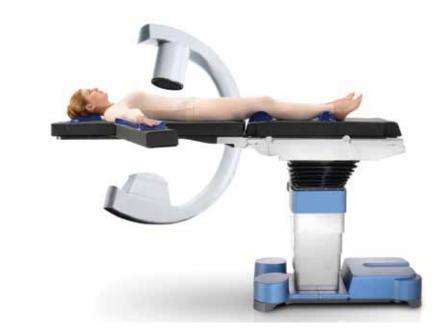







Electro-Hydraulic Operating Table

Flexible Patient Positioning

General Surgery














Lateral position with dual-head plate




More access to patient with dual-joint leg plate

Electro-Hydraulic Operating Table

Technical Specifications

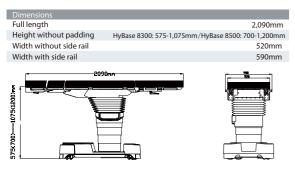

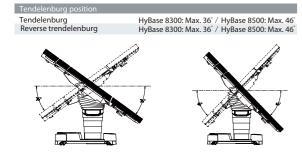

Leg plates up	Max. 80°
Leg plates down (Without sliding)	Max. 90°
Leg plates down (With sliding)	Max. 100°

heliex rosition	
Reflex Position	110°
Trendelenburg	30°
Mounting point up	40°
<u> </u>	
*	



