




PARAMOUNT BED CO., LTD.
2-14-5 Higashisuna Koto-Ku, Tokyo 136-8670, Japan
Tel: +81 3 3648 2961 Fax: +81 3 3648 2420
<http://www.paramount.co.jp/english/>

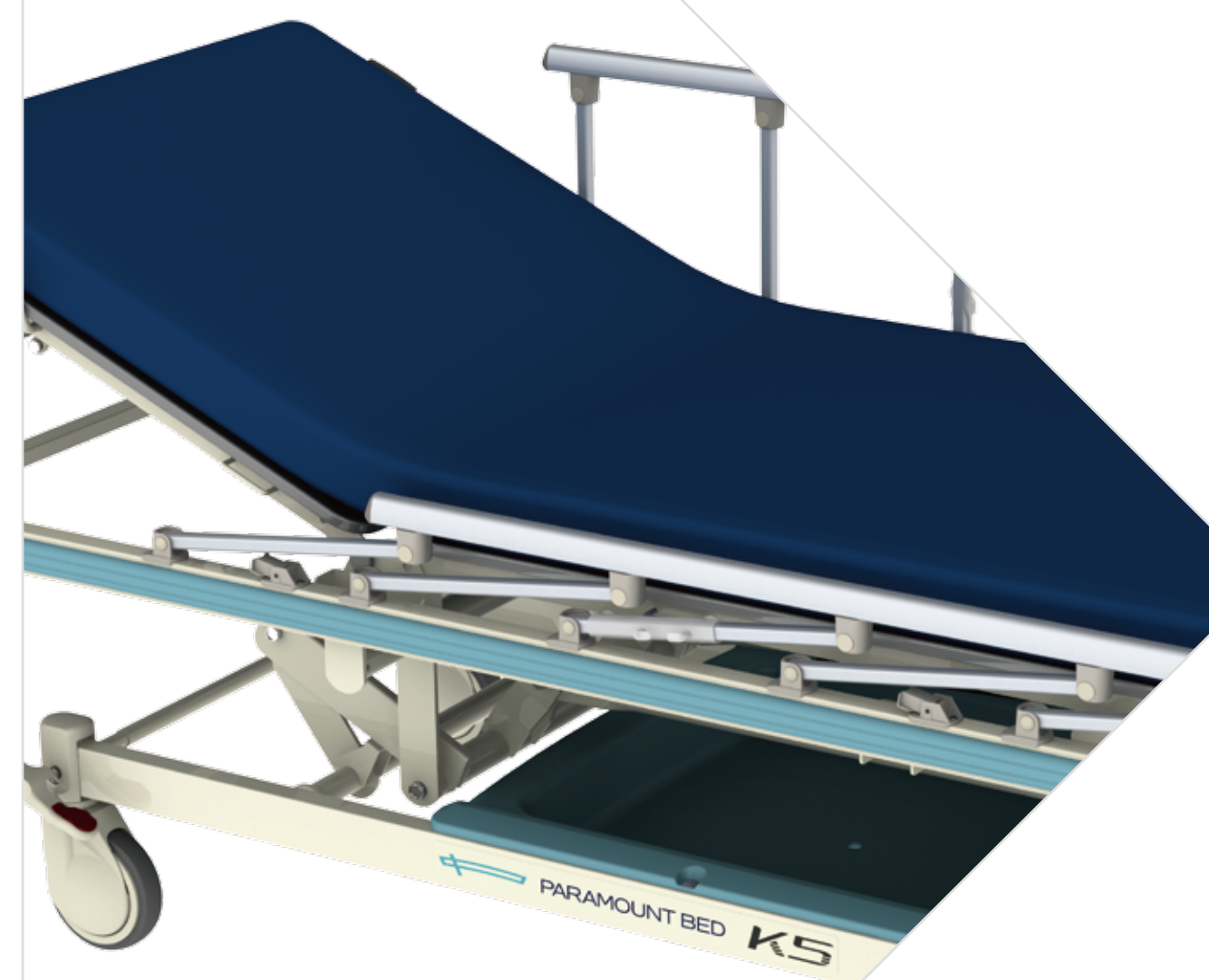


Paramount Bed K5 Series
Stretcher

Designs and specifications are subject to change without notice. All images of the product, its components and accessories may not reflect the final product. The images may reflect optional equipment and are subject to change without notice. All information presented herein is based on current data or design goals available at the time of posting and is subject to change without notice. Some equipment and/or features may be unavailable. Actual results of performance and other specifications may differ or vary with production models and may depend on selected models and/or accessories.

Safety

Side-rail Effective Height

By raising the side-rail effective height, the risk of falling from the stretcher by climbing over the side-rail is reduced.

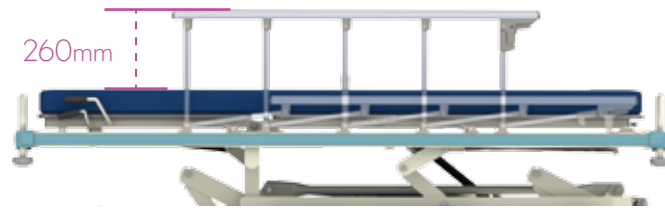

Mattress Fixing Belt

Not only does the mattress fixing belt secure the mattress to the stretcher to prevent slippage, it also secures the patient to the mattress. This is crucial to patient safety when patients are transported on stretchers. In the interests of hygiene, the belt can be removed for washing.


Trendelenburg Reverse / Trendelenburg (PK-5114/5114X only)

In an emergency, trendelenburg / reverse trendelenburg is possible. For safety, the stretcher first stops at the horizontal position when being switched between the trendelenburg and reverse trendelenburg positions.

Work Efficiency

X-ray Imaging

X-ray imaging using the cassette based system can be performed within the range indicated by the shaded area in the diagram.

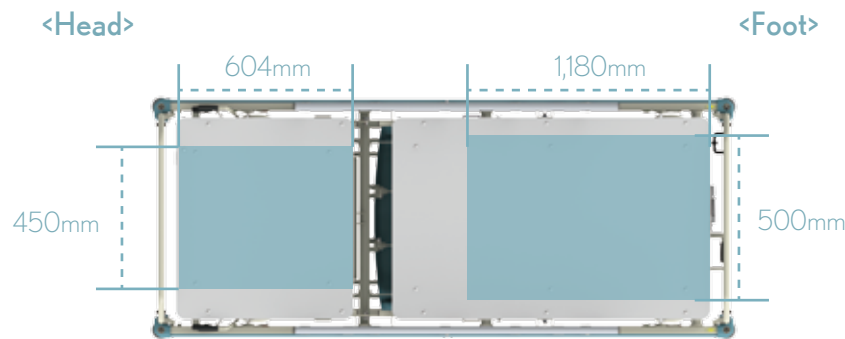

Casters

- 150mm single wheel casters
- Total locking system
 - A preventive measure in situations where medical personnel may miss locking down all casters.
 - Stretcher can be locked down using the caster operation pedal at any of the four casters.

Transport Handle

With transport handles that span the entire width of the stretcher, medical personnel can adjust the position of their arms such that it best suits their shoulder width when pushing the stretcher. This makes stretcher transport easier.





Accessory Hook

To allow miscellaneous items such as a urine bag to be attached, an accessory hook is a standard part with all models. (1 pcs/ each side-rail)
To accommodate different use scenarios, the position of the hooks may be adjusted.

Other Feature

Torque Release Handle

Can be turned easily, which reduces the workload of the caregivers. Also K5 series ' unique torque release handle protects the crank shaft from damage and ensures a long product life.

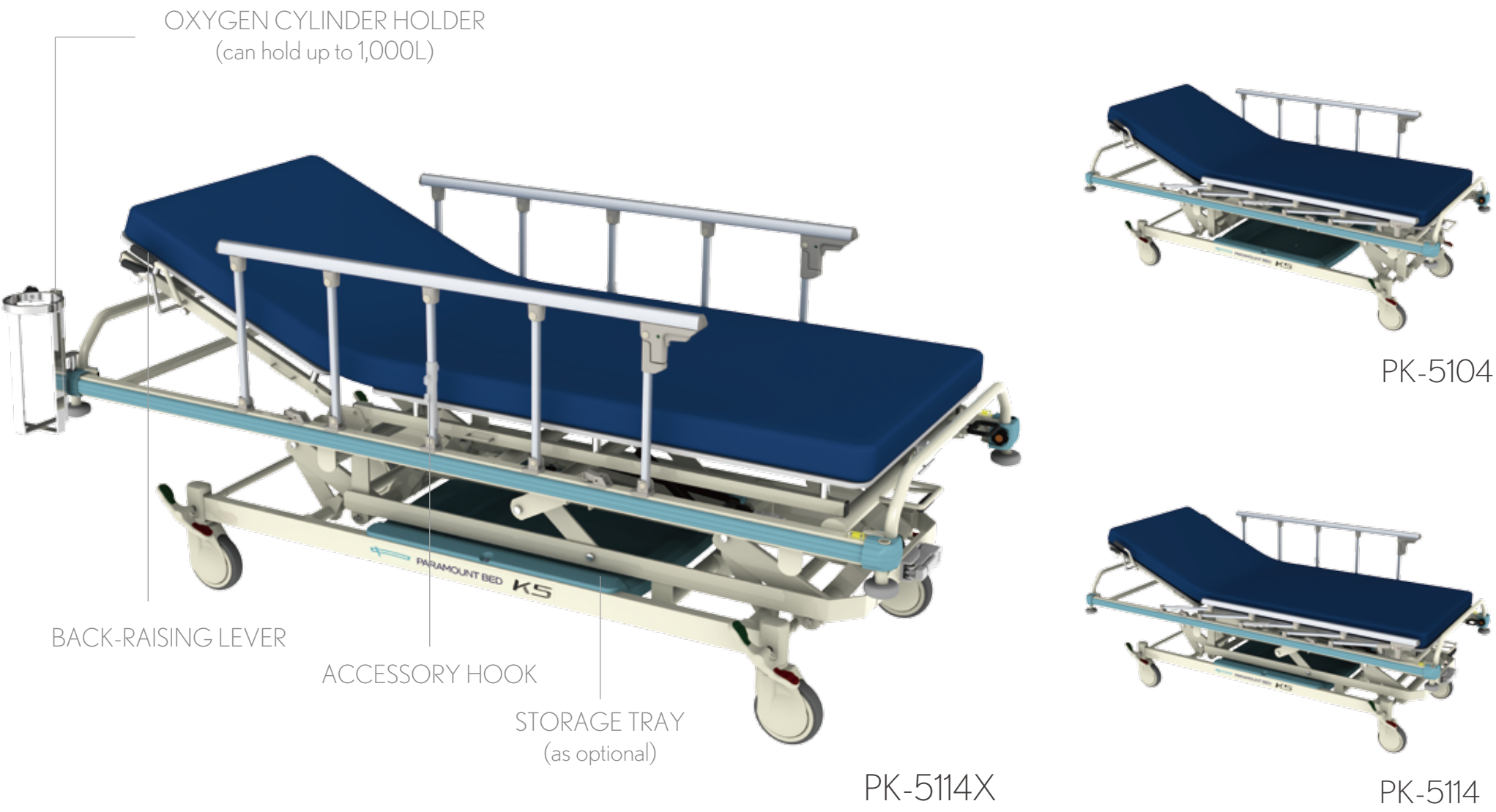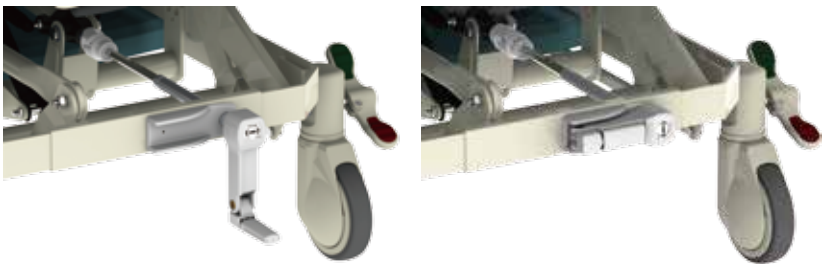

Specifications

Model No.	PK-5104		PK-5114	PK-5114X
Total Length			2,175mm	
Total Width			900mm	
Mattress Base Length			1,985mm	
Mattress Base Width			750mm	
Mattress Base Height*	580 - 900mm			630 - 950mm
Caster	150mm dia. single-wheel casters (total locking system)			
Safe Working Load**	204kgf (2,000N)			
Back Raise Angle	0° to 90°			
Trendelenburg/ Reverse Trendelenburg	—		TR:12°/RTR: -6°	
Accessory Hook			○ (1pcs/ each side-rail)	
Oxygen Cylinder Holder			○	

CE *The mattress base height indicates the height from the floor to the mattress base surface.
**Maximum load with which the bed can be operated. (The total load of user's weight and incidental items including the mattress and options.)

Mattress Features

- MRSA Anti-bacterial Treated
- Flame Retardant
- Water Proofed

Dimension		W750 x L1,985 x Th 80 (mm)
Material	Surface	Vinyl Leather
	Core	Urethane Foam

Additional Accessories

- Foot Board with Writing Table (PK-2100A)
- IV Pole (KC-56A/B)
- Storage Tray (PC-5078)
- Oxygen Cylinder Holder* (PC-5057)

*Can hold up to 500L