

IV Stand KC-508A

A step is provided around the legs that can be depressed with a foot to keep the IV stand from moving or help the caregiver over steps.

* Not stackable.



Optional Accessory

**Oxygen Cylinder Holder
KC-0581 (for KC-508A)**

An auxiliary stay is provided to keep the IV stand from tipping over.



**Infusion Bag Hook
KC-0584 (for KC-508A)**

Attach the hook to the top end of the IV stand support column to enable large - volume infusion bags up to 2,000 cc to be attached at a position closer to the support column than with the previous hook.



* Can also be attached to the former KC-508.

Product name	Model No.	Dimensions (cm)	Product weight	Suitable Table	
				KC-509	KC-508A
IV Stand	KC-509	Legs: 47.9, max. diameter x 139.5 to 207.0, total Height	Approx. 5.5 kg	●	—
	KC-508A	Legs: 47.4, max. diameter x 128.7 to 196.2, total Height	Approx. 4.3 kg	—	●
Oxygen cylinder holder	KC-0591	12.7 (W) x 23.7 (D) x 31.5 (H)	Approx. 1.3 kg	○	—
	KC-0581	12.7 (W) x 23.7 (D) x 41.2 (H)	Approx. 1.5 kg	—	○
Linking arm (for a bed/stretchers)	KC-0582	4.6 (W) x 45.8 (D) x 10.3 (H)	Approx. 500 g	○	○
Linking arm (for a wheel chair)	KC-0583	4.6 (W) x 45.8 (D) x 10.3 (H) (*)	Approx. 665 g	○	○
Infusion bag hook	KC-0584	6.3 (W) x 6.5 (D) x 11.4 (H)	Approx. 60 g	—	○
Urine collection bag holder	KC-058	13 (W) x 13 (D) x 5.3 (H)	Approx. 250 g	○	○

* Not including the width, depth, and height of the holder included with the stand

* Specifications and designs are subject to change without notice.

Local Distributor

PARAMOUNT BED CO., LTD.

Head Office: 14-5 Higashisuna 2-chome, Koto-ku, Tokyo 136-8670, Japan
Tel: +81-3-3648-2961 Fax: +81-3-3648-2420
<http://www.paramount.co.jp/english/>

IV Stand

Easy to use for both the patient and the caregiver.
New Stackable KC-509 Debuts.



Safer and Easier Operation.

Convenient features unique to Paramount Bed.

Infusion Bag Hook

To ensure more stability of the IV stand, large - volume infusion bags up to 2,000 cc can be attached at a position closer to the support column than with the previous hook .



Anti-drop function

The support column comes equipped with an air valve to prevent the support column from dropping suddenly, causing the infusion bag to drop.

Grip can be grasped in a variety of ways

The grip has a unique shape and 12 step angle adjusting function so it can be grasped by the patient in a variety of ways.




Extremely safe height adjustment function

Squeeze the upper/lower stopper grips together to unlock the stopper and let them go to lock it. This method prevents the user from forgetting or failing to tighten the stopper all the way after adjusting the height, which happens often with a fixing knob.



Casters run smoothly and safely

The casters are designed to run smoothly and turn in a small radius. The wheel covers are equipped with protrusions on each side to prevent the wheels from slipping into the gaps of elevator doors. The wheels are also designed so that dust does not build up easily and can be removed for cleaning.

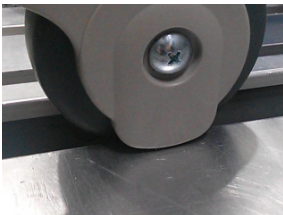


NEW

Stacking saves storage space

IV stands can be stacked together in the horizontal direction. There is a click when they are stacked together correctly so you can verify that they are occupying the minimum space. Multiple IV stands can be stacked together and easily moved around as the casters do not come into contact with the legs.

* KC-509 only

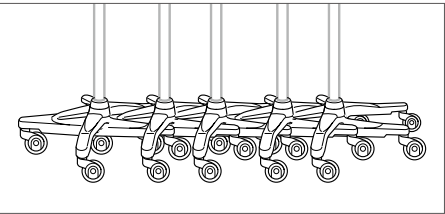

Optional Accessory

Oxygen Cylinder Holder KC-0591 (for KC-509)

An auxiliary stay is provided to keep the IV stand from tipping over.

* IV stands cannot be stacked together when oxygen cylinder holders are attached.



Urine Collection Bag Holder KC-058

Attach the holder to the IV stand support column using a fixing knob to enable the attachment of a urine bag.




Linking Arm

Link the bed, stretcher, or wheelchair to the IV stand with the linking arm to enable them both to be moved together.


The IV stand is linked at a set distance to prevent the patient's catheter getting caught on the IV stand.

KC-0582 (for a bed/stretcher)

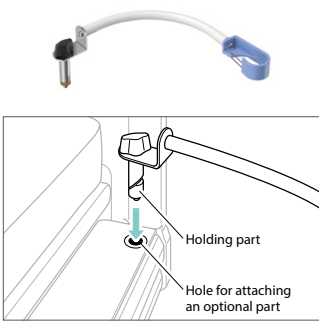

* Option suitable for Paramount beds with option attaching holes, or Paramount Bed Stretcher KK-715/KK-725 series. Secure it to the option attachment hole or IV stand attachment hole.

KC-0583 (for a wheel chair)

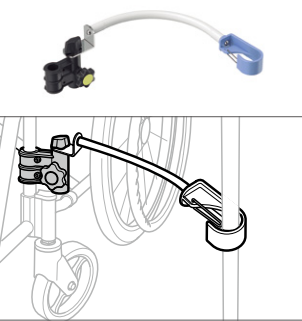

* Option suitable for the Paramount Bed Wheelchair T50/T51/T54 series.
* First attach the holder that comes with the IV stand to the wheelchair and then secure both using the linking arm.