

L'za Series

Multiple Medical Recliner



Suitable for...

- Dialysis**
- Chemotherapy**
- Blood Donor, etc**

- For First Class Care -

***First Class
Comfort***

***First Class
Function***

***First Class
Quality***



PARAMOUNT BED

We promise to provide you “First Class”



POINT 1 For First Class... Comfort


Comfortable Seat Design

The seat is designed very widely (Back Section: 600mm, Leg Section 720mm). Able to freely change the sitting posture, (Ex. Open or cross your legs on the recliner). In addition, by adopting Double Layer Urethane and “S” spring structure which is the same as the advanced seat design for Japanese cars, users can sit with comfort.

Double Layer Urethane (Image)



Wide Seat Design

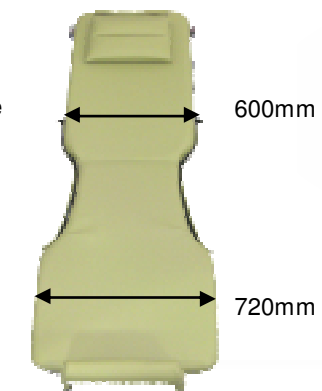

“S” Spring Structure



Comfortable Positions Just Like The Exclusive PARAMOUNT BED

1. Reclining Position

Back and leg sections are adjustable by electric control. Reclining Movement (Auto Contour Movement) is comfortable and luxury for users to rest with relaxation, because we made full use of our advanced technology as the No.1 Japanese hospital bed manufacturer.



2. Fully Reclined Position

“Fully Reclined Position”, which is fully flat from head to foot, is available by L'za. By realizing comfortable sleeping posture, the user can stay with relaxing ease on the recliner, even for long durations. Moreover, by using Reclining button, the position can move from Reclining Position to Fully Reclined Position approx. by 2.3 times* faster.
*Comparison with Back lowering movement.



3. Trendelenburg (Tilting) Position

Trendelenburg (Tilting) Position is available for emergency cases like anemia. L'za is designed for being suitable for various medical uses, Dialysis, Chemotherapy, Blood Donor, etc.



4. Height Adjustment (KA-8260 Series only)

Hi-Lo adjustment is available by electric control. By adjusting Seat height, it is very easy to get on and off. Furthermore, Hi-Lo position is friendly for the injection. Accordingly, safe and efficient care will be provided by using Hi-Low adjustments.



POINT 2 For First Class... Function

Intelligent Hand Held Controller



Universal Design Concept is adopted in consideration for various users. "Easy to hold shape", "Light touch" and "Activating lamps" ensures a safe and friendly use. Integral functions are enabled for caregivers to manage the environment.


Basic Movements (Operated by buttons)

1. Reclining (Auto Contour Raise/Lower)
2. Back Section Raise/Lower
3. Leg Section Raise/Lower
4. Hi-Low Adjustment (KA-8260 Series only)

Integral Functions (Operated by side holes)

- A Reset button
- B 105° Angle Limit and release button
- C Back & Leg Section Raise/Lower speed control button
- D Sound beep and off button
- E Lock-out button for Reclining/Back/Leg/Hi-Low

The picture is 3M.

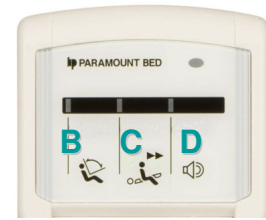


All functions are operated by a pointy tip by pushing holes on side surfaces.

Various Functions Based on Research and Medical Uses

1. Arm Rest

Arm rest can be lifted up in order to get on and off easily, and fixed in the raised position for the easier injection, as well. Moreover, the side adjustable arm rests are suitable for various user's body size.

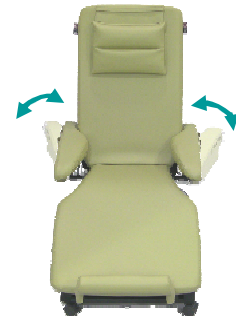

2. Foot Rest

By using the lever on the side of seat, Foot Rest is adjustable. This allows users to place their foot firmly on the foot rest and adjust the sitting position on their own. In case an user gets a cramp in her/his legs, she/he can stretch the muscle of calf by using the Foot Rest.




3. Caster

For the better transportation, 65mm casters (With brakes on each caster) are equipped. Furthermore, 100mm casters are available as an accessory.

REMARK:
When 100mm casters are installed, the height is 51mm higher.

65mm Caster with Brake

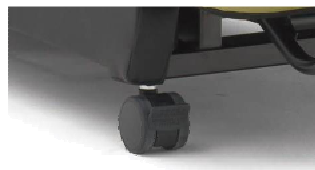


100mm Caster with Brake (Accessory)



4. Pillow

The pillow consists of the pressure relieving urethane. Therefore it is fit and comfortable for user's head. Moreover, the height is adjustable according to the preference.

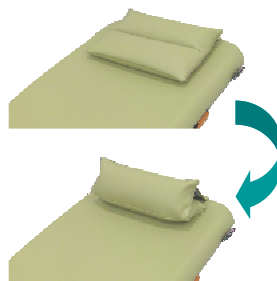

5. E.T. Lever

In emergency cases, E.T. (Emergency Trendelenburg) is available by pulling the orange lever. In addition, in case of power failure, back section can be lowered.



POINT 3 For First Class... Quality

Silent Actuator

It is same as the hospital bed's actuator. Paramount's exclusive Silent Linear Actuator constantly provides the same speed regardless of user weight with less-vibration and silent movement. Electrical safety has been enhanced by complying with all CE regulation.



Seamless Cover

In order to prevent a liquid from sinking into the inside of the cover or accumulating the dust, the joint section of the cover is welded. Moreover, in case the cover gets damaged, the replacement is very easy because the cover is removable.



Accessories

KF-8250 Table



KC-825T TV Arm Holder



REMARK:
If you need TV & TV Arm,
contact our distributor.

KC-825C 100mm Caster



Bed Side Screen (Not CE Certified)

KC-021P



KC-022PE



KC-023BL



KC-56E I.V. Pole KC-825A I.V. Pole Holder



KE-072 Sheets Set

- Washable
- Color: White
- 1 Set is consist of
5 Sheets as below;

Pillow
Main Body
Arm Rests (2 sets)
Foot Rest

KC-826F (For 3M) KC-825F (For 2M) Foot Control Switch

Hands-Free control is
available by Foot Control
Switch. (The picture is
KC-826F.)

KC-508 IV Stand (Not CE Certified)



Specifications

Movements

Back Section
-12° to 75°



Leg Section
-30° to 12°



Height Adjustment
450mm to 760mm
(KA-8260 Series only)




Trendelenburg
0° to 12°



Dimensions

Color variations

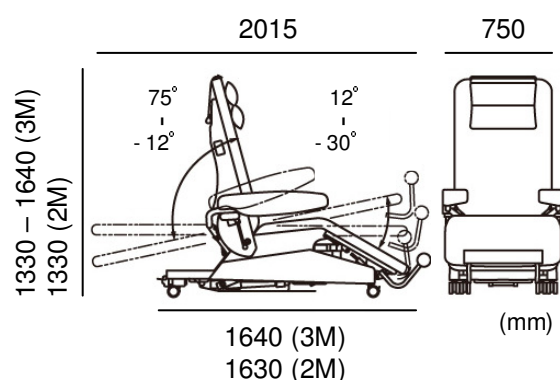

Y: Yellow Green



P: Pink



G: Gray

REMARK: According to preference, please refer to above alphabet code just after the bed model number when ordering. Ex. KA-8260YE (Yellow Green)

Specification Sheet

Description			Multiple Medical Recliner L'za Series	
Model Number			KA-8250YE / KA-8250PE / KA-8250GE	KA-8260YE / KA-8260PE / KA-8260GE
Color Variations			Yellow Green, Pink, Gray	
Electrical Components	Operation Range	Back section	-12° to +75°	
		Leg section	-30° to +12°	
		Height adjustment	-	310mm
Size	Overall length		1630mm	1640mm
	Overall width		750mm	
	Height	Overall Height	1330mm	1330mm - 1640mm
		Seat height	450mm	450mm - 760mm
Caster			Diameter 65mm double wheel caster with stoppers	
Materials	Mainframe		Steel with powder coating	
	Cover		Vinyl Leather	
	Cushion		Urethane Foam	

REMARK: When 100mm casters are used, the height becomes 51mm higher.



PARAMOUNT BED

PARAMOUNT BED CO., LTD.

Head office 14-5 2-Chome Higashisuna Koto-Ku Tokyo 136-8670 JAPAN

Tel: +81-3-3648-2961 Fax: +81-3-3648-2420

<http://www.paramount.co.jp/english/>