

## Lineup

KA-391 without brake type



KA-392 with brake type




Accessory  
KA-039 Basket



Dimensions:  
250 x 169 x 275 (All in mm)

## Specifications

| Product name |                 | Walking Aid                     | Walking Aid (with brake) |
|--------------|-----------------|---------------------------------|--------------------------|
| Product code |                 | KA-391                          | KA-392                   |
| Dimension    | Length          | 732mm                           | 752mm                    |
|              | Width           | 660mm                           | 660mm                    |
|              | Height          | 966 to 1,218mm                  | 966 to 1,218mm           |
|              | Weight          | 12.5 Kg                         | 13.5Kg                   |
|              | Arm rest height | 893 to 1,146mm                  | 893 to 1,146mm           |
| Caster       |                 | 100mm dia. silent caster        |                          |
| Material     | Support pole    | Aluminum (Alumite coating)      |                          |
|              | Base frame      | Steel (Melamine firing coating) |                          |
|              | Arm rest        | Polypropylene                   |                          |
|              | Cushion         | EVA resin                       |                          |
|              | Front bumper    | Polypropylene                   |                          |
|              | Side bumper     | Polypropylene                   |                          |
|              | Rear bumper     | Elastomer resin                 |                          |

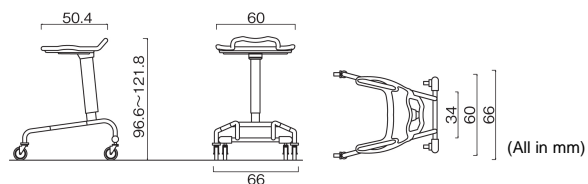

### PARAMOUNT BED CO., LTD.

Head office 14-5 2-Chome Higashisuna Koto-Ku Tokyo 136-8670 JAPAN  
Tel: +81-3-3648-2961 Fax: +81-3-3648-2420  
<http://www.paramount.co.jp/english/>

\*Specifications and designs are subject to change without notice.

2010.6

# Walking Aid

## KA-390 Series



*- The best partner for rehabilitation -*

Easy to  
maneuver

Easy to  
adjust

Easy to  
clean

PARAMOUNT BED

PARAMOUNT BED

## Example of use

KA-390 series are available for each stage, in the recovery process for walking.

Step 1

Step 2

Step 3



*Best partner!*

## Easy to adjust

Every user chooses the preferable height by "Easy to adjust" feature.

*Best partner!*




### Easy height adjustment by high quality gas spring lever

The height of arm rest is easily adjusted according to user's preference, by just pulling the lever. (Between 893mm and 1,146mm freely) The high quality gas spring enhances its usability dramatically.

## Easy to maneuver

"Easy to maneuver" feature makes the user happy to walk.

*Best partner!*

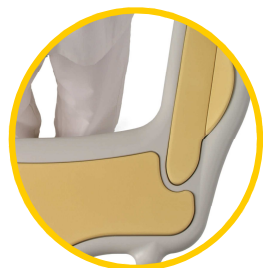

### Grip

The ergonomically designed grip is suitable for various styles of gripping.



### Arm rest

In order to support user's body, the curved-design is very effective because the shape fits to various use, according to user's conditions.



### Elbow rest

Elbow rest makes it easy to operate with peace-of-mind because the rise-edge section supports user's elbow reliably. Also, EVA resin-made cushion reduces the pressure to user's elbow.



### Brake (KA-392 only)

The big-brake lever is very easy to lock/unlock for various users. Moreover, parking brake function is equipped to prevent the walking aid from moving unintentionally.

## Easy to clean

"Easy to clean" feature keeps the walking aid clean for many years.

*Best partner!*



### Detachable Elastomer resin made cushion


The special cushion on arm rest is made of Elastomer resin, which has a chemical proof feature with the elasticity. In addition, the cushions are detachable to clean each part. In case the cushion is needed to change, each part can be exchangeable.

## Other features



### Stacking

The exclusively designed walking aid can save the space to store. Required space for storing is 5 units consumes approx. the spaces of 2 units only.

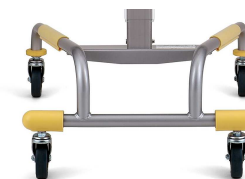

### Corner & side bumper

The corner & side bumpers cover user's feet for safety. Also, the bumpers reduce the damage to wall when the walking aid bumps against.



### 100mm silent caster

Stress-free walking is come true by the 100mm silent casters. Smooth and noiseless walking make the user very comfortable because there is no need to worry about disturbing others.



### Open space of base frame

The front section of base frame is widely opened for easy-walking.