



PT. PARAMOUNT BED INDONESIA  
MM2100 Industrial Town, Block M-1-1 Gandamekar,  
Cikarang Barat, Bekasi, Jawa Barat, 15730, Indonesia  
Tel: 021-8981389 (Direct)  
FAX: 021-8980716  
[www.paramount.co.id](http://www.paramount.co.id)  
[marketing@paramount.co.id](mailto:marketing@paramount.co.id)



K-HOUSE (JAKARTA SHOWROOM)  
Jl. Prof. Dr. Satrio No.22 Jakarta Selatan Tel: 021-5200716  
Opening Hours: Mon-Fri 9a.m. - 5p.m.



Designs and specifications are subject to change without notice. All images of the product, its components and accessories may not reflect the final product. The images may reflect optional equipment and are subject to change without notice. All information presented herein is based on current data or design goals available at the time of posting and is subject to change without notice. Some equipment and/or features may be unavailable. Actual results of performance and other specifications may differ or vary with production models and may depend on selected models and/or accessories.

L183-3273-EN01



PA-8200 Series  
Multiple Medical Recliner

designed for  
dialysis  
treatment

# Product Features Highlight

## Handheld Controller

Handheld controller features a slim width ergonomically sculpted grip, making it easier to clean.



## Emergency Back Lowering Lever

This function is to lower the backrest so it can be used rapidly in case of emergency, such as when patient is fainting due to anaemia or a similar condition or a power failure.



## Seamless Seat Cover

In the interest of hygiene, the seat cover is made with a seamless material that is impermeable.

To clean soiled seat covers, simply wipe away the contaminating liquid.




## Armrest Movement Flexibility

Armrest can be adjusted to move up and down by using plunger to make approach and getting on or off the recliner easily.



## USB Port\*

USB Port are located on the right side of base frame. It can be used for charge the battery of mobile phone, PC, etc.

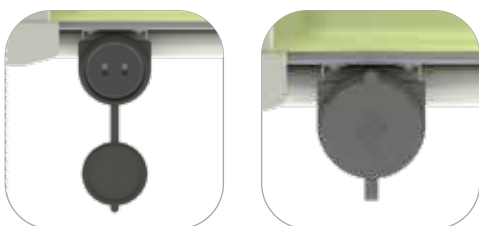

\*Only for PA-8210

## Position



Back Raise & Lower Leg Position

Flat Position





Product Code		PA-8200	PA-8210
Motors		2 motors	
Total Length		1,480mm (When flat 1,900 mm)	
Total Width		980mm	
Total Height (*1)		1,370mm	
Seat Width		550mm - 720mm	
Seat Height		570 mm	
Casters		100mm dia with stopper	
Power Cable Length		4,000 mm	
Material	Mattress Base	Plywood with melamine	
	Main Frame	Made of Steel, Powder Coating	
	Arm Rest	Vinyl Leather, Urethane Foam	
User' s Maximum Weight		138 kg	
Safe Working Load		204kgf ( 2,000 N )	
Back Raise Angle		0° to 75°	
Leg Raise Angle		-30° to 0°	
USB Port		—	○

\*1) Dimension when Back Raise Angle Maximum

## Additional Accessories

- Roll Sheets Holder ( PC-5104 )
- IV Pole ( KC-56A/B )
- IV Pole Attachment Hole ( PC-5103 )
- Pillow ( PE-8201 )
- Foot Rest ( PC-5102 )



### Adjustable Pillow Height

To Make easier for the patient to adjust the position of the pillow with one hand, a free-stop pillow with a weight for balancing is used, and the pillow is attached by a belt carrying the weight.



### Easily Adjustable Foot Rest

There is two holes on the frame side that can be used to adjust footrest position according to height patient by using knob bolt.

