



Neuron MAX™

6 F Long Sheath

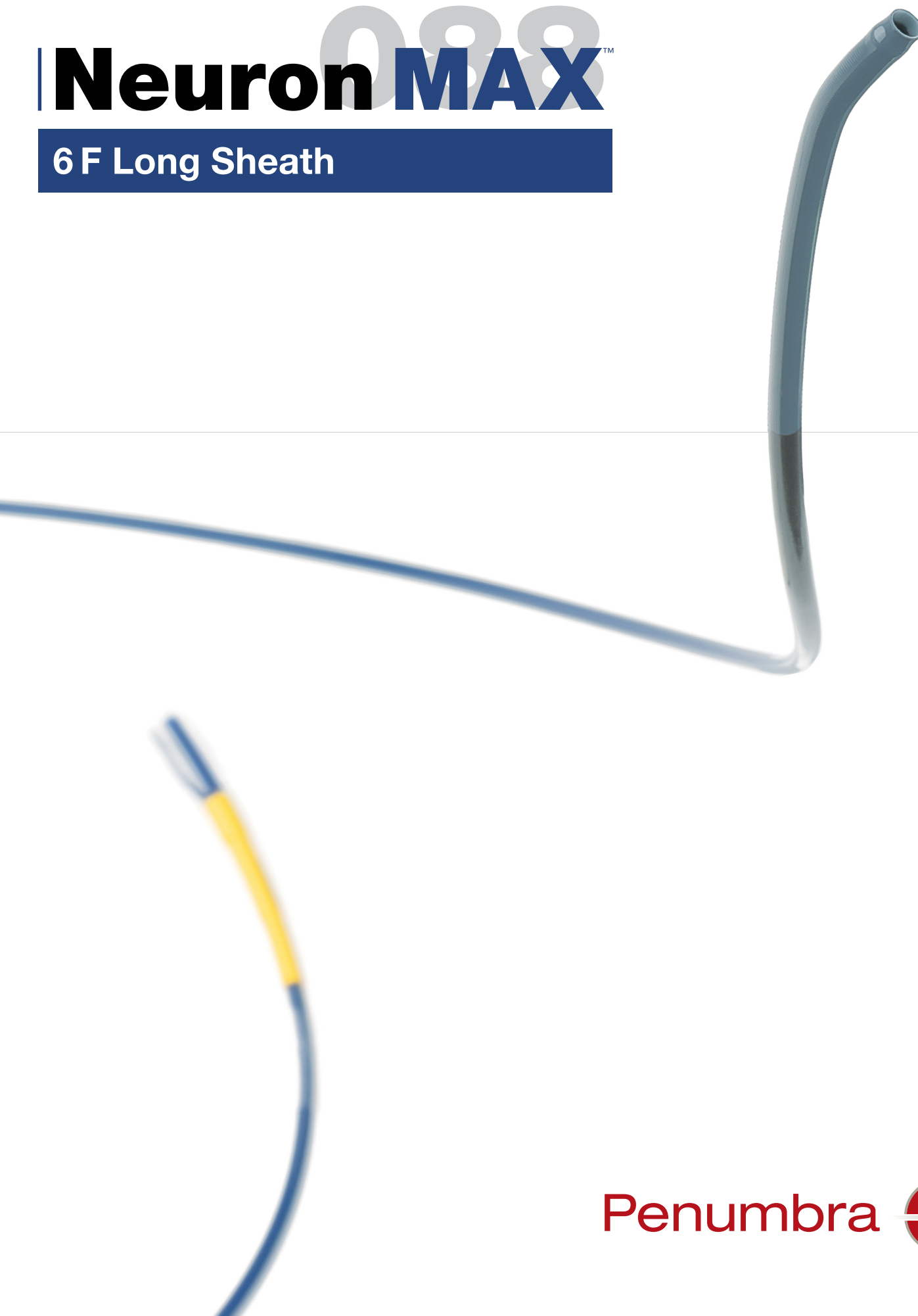

Neuron MAX 088 Advantages

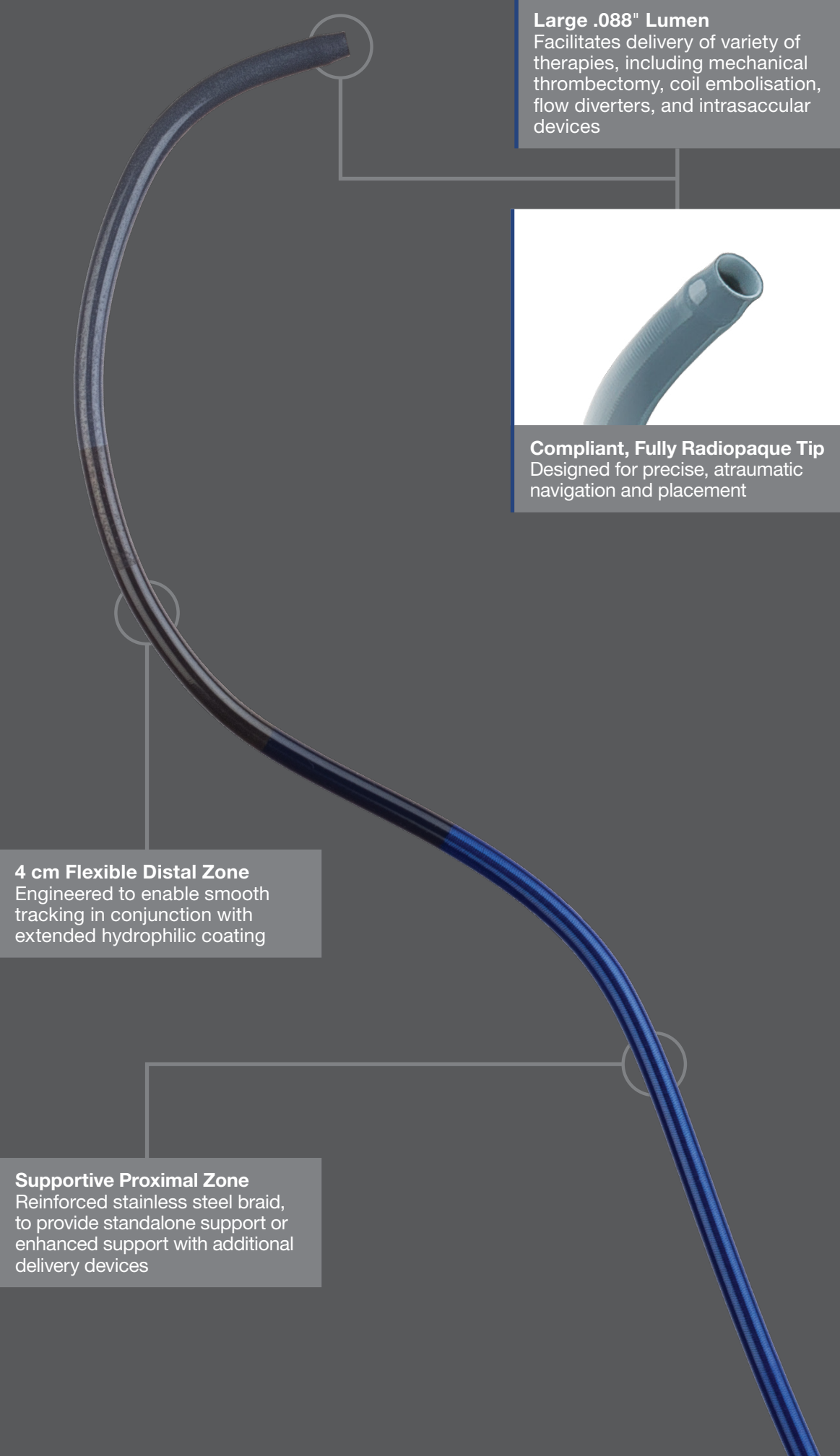

Atraumatic, Radiopaque Tip

Neuron MAX 088
Full Distal Shaft Radiopacity + Marker Band

Stryker® AXS Infinity LS® PLUS™
Radiolucent Tip + Marker Band


Stryker AXS Infinity LS
Radiolucent Tip + Marker Band

Neuron MAX 088 vs. Stryker AXS Infinity LS PLUS


Neuron MAX 088 vs. Stryker AXS Infinity LS

Devices shown at same scale. Photographs taken by and on file at Penumbra, Inc. Tests performed and data on file at Penumbra, Inc. Bench test results may not be indicative of clinical performance.

Distal Flexibility

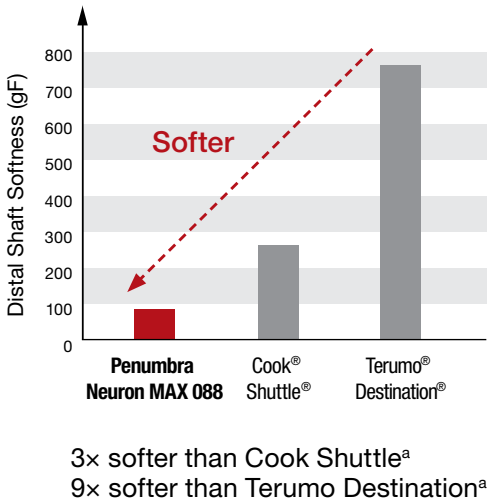

Optimal Proximal Support

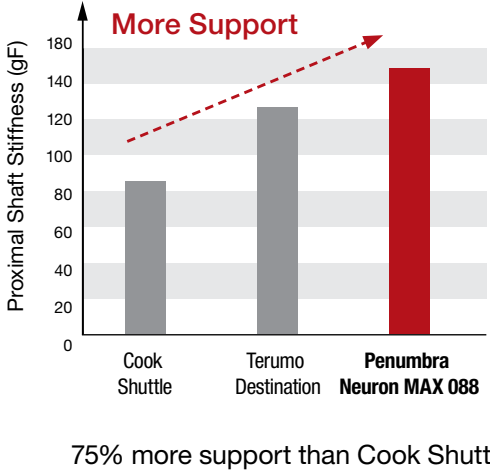

Typical Neuron MAX 088 Delivery

Designed to be positioned at high cervical ICA (or more distal) for maximum stability

1

6 F Select™ Catheter
125 cm, 6.5 F max OD

.035 / .038" Wire

Seamless transition zone designed for atraumatic advancement

Neuron MAX 088
90 cm

2

Typical **Neuron MAX 088** Placement

Rendering for illustrative purposes only. Individual results may vary depending on patient-specific attributes and other factors.


a. Tests performed and data on file at Penumbra, Inc. Bench test results may not be indicative of clinical performance.

Designed for Rapid Primary Access


Neuron MAX™ 088 is Compatible with Penumbra 6 F Select™ Catheter to enable Rapid Primary Access

- Facilitates rapid primary access and atraumatic placement into desired vessel without over-the-wire exchange
- .040" lumen accommodates diagnostic angiograms
- Stainless steel braided proximal shaft engineered for support and torque response, with a soft polymer distal shaft to enable atraumatic vessel selection

6 F Select Catheter | 105 / 125 cm

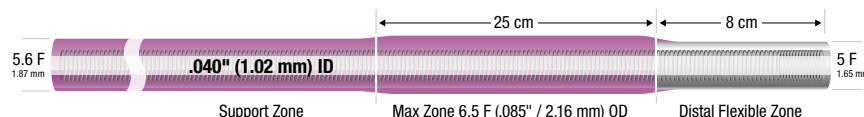

Typical Penumbra System™ Stroke Setup Using Neuron MAX 088

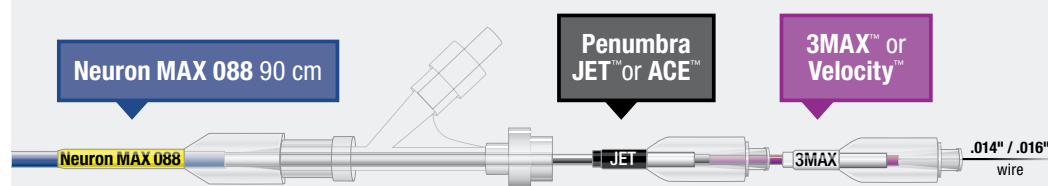

Penumbra ENGINE™
with Penumbra
ENGINE Canister

Reperfusion
Catheter

Hi-Flow Aspiration Tubing

Ordering Information

Catalog Number	Description	Working Length (cm)	Distal Flexible Zone (cm)	OD (Proximal / Distal) (F)*	ID (in.)	Wire Compatibility (in.)
Neuron MAX 088 6 F Long Sheath (Crosscut Valve, RHV, and Dilator Included)						
PNML6F088804	Neuron MAX 088 6 F Long Sheath, 80/4 Straight	80	4	8 / 8	.088	.035 / .038
PNML6F088804M	Neuron MAX 088 6 F Long Sheath, 80/4 MP	80	4	8 / 8	.088	.035 / .038
PNML6F088904	Neuron MAX 088 6 F Long Sheath, 90/4 Straight	90	4	8 / 8	.088	.035 / .038
PNML6F088904M	Neuron MAX 088 6 F Long Sheath, 90/4 MP	90	4	8 / 8	.088	.035 / .038
PNML6F0881004	Neuron MAX 088 6 F Long Sheath, 100/4 Straight	100	4	8 / 8	.088	.035 / .038
PNML6F0881004M	Neuron MAX 088 6 F Long Sheath, 100/4 MP	100	4	8 / 8	.088	.035 / .038
6 F Select Catheters						
PNS6F105BER	6 F Select Catheter, 105 BER	105	9	5.6 / 6 / 5	.040	.035 / .038
PNS6F105H1	6 F Select Catheter, 105 H1	105	9	5.6 / 6 / 5	.040	.035 / .038
PNS6F125BER	6 F Select Catheter, 125 BER	125	9	5.6 / 6 / 5	.040	.035 / .038
PNS6F125SIM	6 F Select Catheter, 125 SIM	125	9	5.6 / 6 / 5	.040	.035 / .038
PNS6F125H1	6 F Select Catheter, 125 H1	125	9	5.6 / 6 / 5	.040	.035 / .038
PNS6F125SIMV	6 F Select Catheter, 125 SIM-V	125	9	5.6 / 6 / 5	.040	.035 / .038

*A mid outer diameter is only listed if applicable to device.

Prior to use, please refer to the Instructions for Use for complete product indications, contraindications, warnings, precautions, potential adverse events, and detailed instructions for use.

NEURON MAX System – Intended Use

The NEURON MAX System is intended for the introduction of interventional devices into the peripheral, coronary, and neuro vasculature. **Potential Adverse Events** Possible complications include, but are not limited to, the following: acute occlusion; air embolism; death; distal embolization; emboli; false aneurysm formation; hematoma or hemorrhage at puncture site; infection; intracranial hemorrhage; ischemia; neurological deficits including stroke; vessel spasm, thrombosis, dissection, or perforation.

PENUMBRA SYSTEM™ – Intended Use

The PENUMBRA SYSTEM is intended for use in the revascularization of patients with acute ischemic stroke secondary to intracranial large vessel occlusive disease using continuous aspiration. **Potential Adverse Events** Possible complications include, but are not limited to, the following: allergic reaction and anaphylaxis from contrast media; acute occlusion; air embolism; arteriovenous fistula; death; device malfunction; distal embolization; emboli; false aneurysm formation; hematoma or hemorrhage at access site; inability to completely remove thrombus; infection; intracranial hemorrhage; ischemia; kidney damage from contrast media; neurological deficits including stroke; vessel spasm, thrombosis, dissection, or perforation; radiation exposure that may lead to cataracts, skin reddening, burns, alopecia, or neoplasia from x-ray exposure.

PENUMBRA SYSTEM – Intended Use

The PENUMBRA SYSTEM is intended for use in the

revascularization of patients with acute ischemic stroke secondary to intracranial large vessel occlusive disease using continuous aspiration. **Potential Adverse Events** Possible complications include, but are not limited to, the following: allergic reaction and anaphylaxis from contrast media; acute occlusion; air embolism; arteriovenous fistula; death; device malfunction; distal embolization; emboli; false aneurysm formation; hematoma or hemorrhage at access site; inability to completely remove thrombus; infection; intracranial hemorrhage; ischemia; kidney damage from contrast media; neurological deficits including stroke; vessel spasm, thrombosis, dissection, or perforation; radiation exposure that may lead to cataracts, skin reddening or burns from x-ray exposure.

PENUMBRA ENGINE – Intended Use

The PENUMBRA ENGINE is intended as a vacuum source for Penumbra Aspiration Systems.

Penumbra Delivery Microcatheters – Intended Use

The Penumbra Delivery Microcatheters are intended to assist in the delivery of diagnostic agents, such as contrast media, and therapeutic agents, such as occlusion coils to the peripheral and neuro vasculature. **Potential Adverse Events** Possible complications include, but are not limited to, the following: acute occlusion; hematoma or hemorrhage at access site; death; intracranial hemorrhage; hemorrhage; infection (at access site); distal embolization; ischemia (cardiac and/or cerebral); embolus (air, foreign body, thrombus, plaque); aneurysm perforation; false aneurysm formation; neurological deficits including stroke; vessel spasm, thrombosis, dissection, perforation or rupture; air embolism; emboli.

Tip Shapes

Neuron MAX 088

Select Catheter



Penumbra, Inc. USA
One Penumbra Place
Alameda, CA 94502
USA
1.888.272.4606
T 1.510.748.3200
F 1.510.748.3232
order@penumbrainc.com
info@penumbrainc.com

Penumbra Europe GmbH
Am Borsigturm 44
13507 Berlin
Germany
T +49 30 2005 676-0
F +49 30 2005 676-10
de-order@penumbrainc.com
de-info@penumbrainc.com

Penumbra Neuro
Australia Pty Ltd
Suite 3, Level 5, 1 Oxford Street
Darlinghurst NSW 2010
Australia
T +61-1300 817 025
F +61-1300 817 026
order.anz@penumbrainc.com

Penumbra Latin America
Distribuidora de Equipamentos
e Produtos Médicos Ltda
Avenida Brigadeiro Luís Antônio
3421 cj 201 CEP 01401-001
São Paulo, Brazil
T 5511.2883.5825
order.la@penumbrainc.com

Product availability varies by country. Rendering for illustrative purposes only. Please contact your local Penumbra representative for more information.

Copyright ©2021 Penumbra, Inc. All rights reserved. The Penumbra P logos, Neuron, MAX, Neuron MAX, Select, Penumbra System, Penumbra JET, ACE, Velocity, and Penumbra ENGINE are registered trademarks or trademarks of Penumbra, Inc. in the USA and other countries. All other trademarks are the property of their respective owners. 12356, Rev. A 03/21 OUS

