



Aeon7200A

Anesthesia Machine

CE 0123

AEOMED
Reliable Quality Thoughtful Service



Classic & Reliable

- A classic anesthesia machine approved by market through the world
- CE certified to provide safe and efficient anesthesia solution
- Integrated sockets for other medical equipment
- Optional AGSS provides more protection for OR environment.
- Optional EtCO₂ and multiple anesthetic gas monitoring (AGM).

Smoothing operation experience

- Easier operation and observation with 8.4" LCD color touch screen and waveforms & parameters in different colors.
- Supporting two vaporizers with interlock to ensure the safety while giving more flexible choices.
- Double-tube mechanical flowmeter for each gas is suitable for various requirements such as Low-flow.
- APL valve with fast release function improves the efficiency of induction and resuscitation.
- One-Key Switch from manual to mechanical ventilation.






Reliable breathing system

- Integrated sensors, non-consumable design
- Large capacity CO₂ absorber which is easy to disassemble
- Effectively reduce the capacity of the system and increase the concentration of anesthetic
- Low resistance inspiration and expiration valves make it more sensitive to support pediatric ventilation

Excellent ventilation management

Comprehensive ventilation modes ensure flexibility in the different anesthetic phases

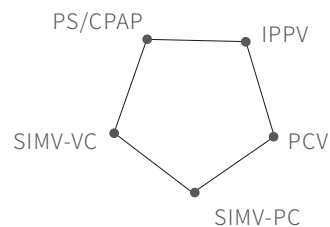

The system monitors the EtCO₂ values and adjusts the parameter to maintain the set EtCO₂ target values.



Spirometry loops and waveforms make it more intuitive on showing the ventilation effectiveness.



Lung recruitment maneuvers include sustained inflation and PEEP increment.

Technical Specifications

Physical & Environmental Specifications

• Dimensions(H×W×D):	1356×845×620mm
• Weight:	110kg (Basic unit)
• Storage temperature:	-20~55°C
• Working temperature:	10~40°C
• Storage humidity:	≤93%

Electrical Specifications

• Mains:	AC 100~240V, 50/60Hz
• Battery:	DC 24V, minimum 120min
• Mains outlet:	3

Pneumatic Specifications

• Gas supply:	0.28~0.6MPa
• Flowmeters(basic config):	O ₂ : 0~1L/min; 1~10L/min N ₂ O: 0~1 L/min; 1~10L/min
• Flowmeters(Full config):	O ₂ : 0~1L/min; 1~10L/min N ₂ O: 0~1 L/min; 1~12L/min Air: 0~1 L/min; 1~15L/min

• Gas system:	O ₂ supply failure alarm O ₂ supply failure protection O ₂ flush: 25~75L/min
• Working mode:	Closed, Semi-closed, Semi-open
• Driven mode:	Pneumatically driven and electronically controlled
• Safety valve:	≤10kPa
• Operating mode:	Man. / Vent.
• Cylinder yokes:	Optional

Ventilator Specifications

• Patient type:	Adult, pediatric
• Ventilation modes(basic config):	IPPV, PCV, SIMV-VC, Manual, PS/CPAP (with apnea backup), Standby

• Ventilation modes(full config):

IPPV, PCV, SIMV-VC, SIMV-PC, Manual, PS/CPAP (with apnea backup), Standby
--

• Setting:

Navigator knob & touchscreen

• Tidal volume:

20~1500ml

• Ventilation frequency:

2~100bpm

• I:E ratio:

4:1~1:8

• Electronic PEEP:

OFF, 3~30cmH ₂ O

• Inspiratory plateau:

OFF, 5~60%

• Flow trigger:

1~15L/min

• Pressure range:

5~70cmH ₂ O

• Pressure support:

3~50cmH ₂ O

• Ventilation Monitoring: Pressure: Ppeak, Pplat, Pmean, PEEP

Volumes:

MV, Vt

Respiratory frequency

Compliance, etc.

• Gas Monitoring: O₂, CO₂, Agent (All optional)

• Graph display(basic config): P-t, F-t, V-t waveforms

• Graph display(full config): P-t, F-t, V-t waveforms,

P-V and F-V loops

• Alarm: High/Low Airway Pressure

High/Low Minute Volume

High/Low FiO₂ (when FiO₂ function is enabled)

Power failure, O₂ Supply failure, Apnea, etc.

Vaporizer

• Support 2 vaporizers (Selectatec® with interlock)

• Optional additional 1 vaporizer holder

• Agent type: Halothane, Enflurane, Isoflurane, Sevoflurane

• Filling type: Pour-fill, Key-fill, Quik-fill®

Remark: Above configurations include standard and option. Please check price with your Aeonmed sales representative.

AEONMED
Reliable Quality Thoughtful Service

The information contained in this leaflet is correct at the date of publication. The policy of Beijing Aeonmed Co., Ltd. is one of the continued improvement to its products. Because of this policy, Aeonmed reserves the right to make any change, which may affect the information in this leaflet without giving prior notices.

HQ: Building 9, No.26 Outer Ring West Road, Fengtai District, Beijing 100070, China
Science Park: No. 10, Chaobai Street, Yanjiao Development Zone, Sanhe City, Hebei Province
065201, China
TEL: +86-10-5841 1198
Http: //www.aeonmed.com/en
FAX: +86-10-8368 1616-8130
E-mail: marketing@aeonmed.com

AAM7200A-2022.03